

RXP 210: Air-volume alarm unit

How energy efficiency is improved

Areas of application

Monitoring laboratory fume cupboard as per EN14175 in combination with an air-volume controller, type RLP100 F123, and the RXP 212 fume cupboard operating unit. Zone 1 potentially explosive areas.

Features

- Comparison of sash opening setpoint with actual value of air volume
- Feeding setpoint difference to alarm unit
- Acoustic alarm in the event that the actual value falls below or exceeds the specified value
- Thermoplastic housing suitable for wall or top-hat rail mounting (rail EN 60715)
- Compressed-air connections with Rp 1/8" female thread
- Low-pressure connections in form of stepped nipples for flexible plastic hose (internal Ø 4 and 6 mm)
- ATEX certification for use in Zone 1 potentially explosive areas
- Conformity tested as per EN 13463-1 and EN 1127-1 (Ex II 2 G T6)
- Complies with directive 97/23/EC Art. 3.3 on pressure equipment

Technical description

- Supply pressure 1.3 bar ± 0.1
- Response sensitivity 1%

Type	Description	Air capacity ³	Air consumption ¹⁾	Weight kg
RXP 210 F001	for monitoring the air flow	180 l _n /h	40 l _n /h	0,2
Supply pressure ²⁾	1,3 bar ± 0,1	Response sensitivity		1%
Input pressures		Perm. ambient temp.		0...55 °C
connector 3	0...1,3 bar	Degree of protection		IP 20
connector 6	0,1...0,6 bar	Connection diagram		A09890
Output pressures		Dimension drawing		M297240
connector 2	0...1,3 bar	Fitting instructions		MV 505381
connector +	approx. 1,0 bar			

Accessories

- 0296936 000*** Bracket for rail EN 60715, 35 × 7,5 and 35 × 15
- 0297354 000*** Short screw-in connector R 1/8 for soft plastic tubing of internal dia. 4 mm, four pieces required
- 0297797 001*** Manometer with indication of deviation in air volume

*) Dimension drawing or wiring diagram are available under the same number

- 1) Extra air consumption when whistle is activated: 180 l_n/h
- 2) For regulations on the quality of the air supply, especially at low ambient temperatures, see section 60
- 3) Air output at connector +

Operation

Together with the RLP 100 F123 controller, the alarm unit monitors the air volume in the fume cupboard. The amount that the fume cupboard's sash is open is fed as a setpoint to the relevant VAV controller and compared to the actual value of the air volume. Any deviation is indicated visually in accordance with DIN 12924 (RXP 212 operating unit). In the event of a malfunction, an acoustic signal is emitted, e.g. whenever the actual value of the air volume is either smaller (min.) by > 8% $\dot{V}$ or larger (max.) by 10% $\dot{V}$ than the setpoint (factory setting). The 'min.' adjuster should not be set to a difference of > 8% $\dot{V}$. The whistle sounds after a delay of four to twelve seconds. It can be stopped by closing connector 6 (pressing the 'mute' button on the RXP 212 operating unit for about two seconds). When actual value again equals setpoint, the acoustic alarm is again ready for operation.


Technical information

Technical manual VAV 7 000621 003




67.531/2 RXP 210

Connection diagram


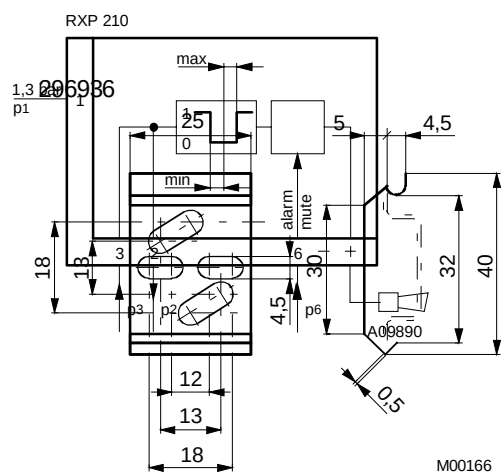


M00166

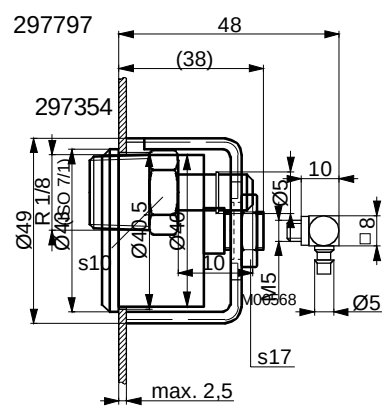


M04020

Dimension drawing

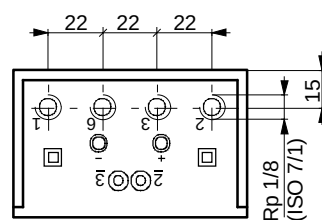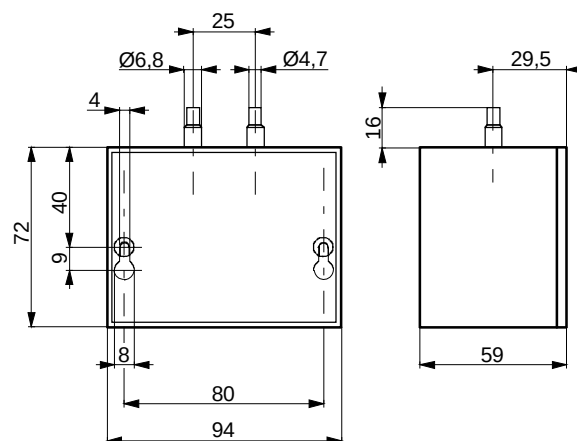

M297240c

Alarm whistle (included in delivery)

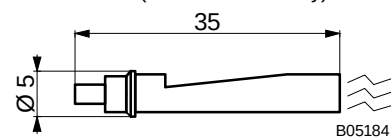

B05184

